

EMMA™ Capnograph

Portable Real-time Capnography



- ❖ CONFIRMS ENDOTRACHEAL TUBE PLACEMENT
- ❖ ETCO2 IN FIRST FEW BREATHS

- > **Immediate results**— minimal warm-up time, with end-tidal carbon dioxide (EtCO₂) and respiration rate (RR) measurements and real-time EtCO₂ waveform displayed in 15 seconds
- > **Small, portable capnograph**— light weight design fits in the palm of your hand for mobility and convenience during short-term EtCO₂ monitoring of adult, pediatric, and infants
- > **Designed to fit easily onto a breathing circuit** — flexible use at multiple points of care including pre-hospital transport, emergency room, operating room, intensive care unit
- > **Rugged design** — for reliable operation in challenging environments
- > **Easy to maintain** — no routine calibration required

Features

- > **Clear & continuous**
Capnogram of carbon dioxide values
- > **Simple, easy-to-use interface**
For quick set-up and one-touch programming
- > **Audible and visual alarm system**
For no adapter, clogged adapter no breath (Apnea), low battery and adjustable high and low EtCO₂ alarm
- > **Battery life**
For up to 10 hours of normal use with two standard AAA lithium batteries

Alarm Limits

Alarm limits and visual alarm status indicator, silence active alarms for two minutes

Capnogram

14.4 second sweep of CO₂ values

Power Button

Warm-up time to full accuracy in 15 seconds

Adapters

1. Light weight
2. Negligible dead space

End-tidal Carbon Dioxide

Quantitative EtCO₂ is updated every breath

Respiration Rate

RR is displayed after two breaths and updated every breath

Airway Adapter

Available in adult/pediatric and infant sizes

Kit




- > **EMMA Kit***
EMMA (mmHg) PN 3639
Net weight 60 grams

Accessories



- > **EMMA Airway Adapter**
Infant
Box of 10
PN 17449

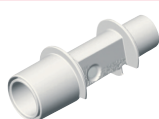

- > **EMMA Airway Adapter**
Adult/Pediatric
Box of 25
PN 17448

Contact Us

+91 70367 70007

+91 96762 27222 / +91 93916 89098

For more information visit our website : www.srisan.in

Flat.No.211, Aparna Green, Opp. Golf View Apartments, Nanakramguda, Gachibowli, Serilingampally, Hyderabad, Telangana, India 500 032.